

UNC STREET DRUG ANALYSIS LAB

How to use the kits:

Our easy-to-use kits allow you to collect a drug sample using a pinhead of powder, a sliver of a pill, or a used cotton. Our machines are super sensitive, so we can also run analyses off of residue sticking to the inside of baggies.

Kits Include:

- Plastic Bag
- Parafilm
- Drapes
- Pencil
- Swab
- Spatula
- 10mg scoop
- 1.5mL of methyl cyanide (acetonitrile)

*This is an organic solvent that renders the sample unusable. Allows us to offer the service through the mail.

After Testing


You return the kit to us, we run it in the lab via GC-MS, and you can check the results here on our website: streetsafe.supply/results

It takes a week (or less) to run the samples once we get them in the mail. Rush service available during outbreaks.

- 1 **Unfold cloth, lay out supplies.**
- 2 **Unpeel tape, unscrew vial.**
- 3 **Wear gloves to prevent contamination.**



video

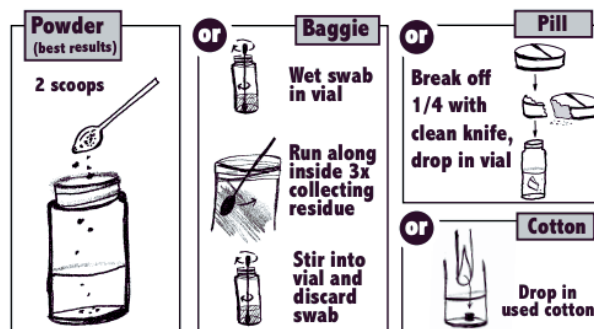

Syringe, foil, and pipe residue are not ideal for analysis.

- 4 **Close vial tightly. Seal using wax tape.**

Remove tape from white backing.



- 5 **Complete back of this card. Card and vial go back in bag.**

- 6 **Give QR code to donor for results.**

Enter Sample ID for results at <https://streetsafe.supply>

Our lab/program is a national service for public health organizations. As a non-profit, we can provide at-cost services to health departments, clinics, and universities. We operate on a per-sample price, offering steep discounts for harm reduction programs as well as free services for drug user unions. Law Enforcement is prohibited from using our service.

www.streetsafe.supply

Learn more at
opioiddata.org